



Cell culture media
For better lives

CELLiST™ CHO Perfusion Media_P1

Shortcut to
Official Channel



CELLiST™ CHO Perfusion Media_P1

I Overview

Ajinomoto Group collaborated with Chugai Pharmaceutical to improve productivity and culture longevity in perfusion culture. Through our joint development project, we optimized both the culture process and the medium composition. By focusing on refining the medium and process, we achieved significant advancements, resulting in higher productivity and longer culture longevity. Our proprietary technologies, such as refined iron sources and exclusive Cysteine technology, position us as leaders in enhancing perfusion culture productivity.

I Key Features of CELLiST™ CHO Perfusion Media_P1

- Easy to prepare and use, easily dissolves at neutral pH.
- Facilitates sustained high cell density and viability throughout the perfusion process.
- Maintains appropriate cell growth rate, reducing requirement for excessive cell bleeding.
- Optimized for suppressing cell death, reducing membrane fouling.

I Properties

- Chemically defined, protein-free medium without any animal-derived components.
- Suitable for all CHO cell lines and their derivatives.
- High performance in both cell growth and protein production.
- Test samples as well as bulk size orders are available.
- Flexible application for easily replacing any existing media platform.
- Manufactured in a GMP-compliant factory.




I Background

The advent of perfusion culture has addressed many challenges faced in batch or fed-batch culture, such as short cultivation periods and purification issues due to high cell density. As a result, biopharmaceutical companies are increasingly adopting perfusion culture. CELLiST™ has developed high-performance perfusion medium for CHO cells, designed for high productivity beyond fed-batch processes. Perfusion media_P1 maintains high cell viability and productivity with a low bleeding rate over extended periods, thanks to exclusive Cysteine technology. This cost-effective solution is ideal for long-term perfusion processes, leveraging Ajinomoto Group's media science expertise.

I Specifications

Format	Powder
Concentration	37.4 g/L
Includes	Poloxamer
Storage Condition	Store in refrigerated (2-8C), dark conditions
Shelf Life	12 months
Item Description	Powder, chemically-defined and animal origin free medium.

I Media Performance

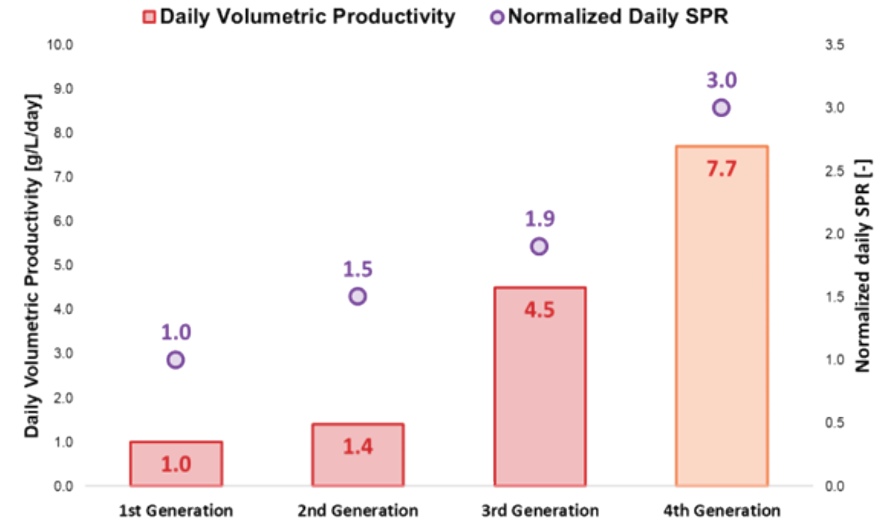

- Increasing poloxamer concentration enhanced cell-specific productivity and viability by protecting cells from shear stress.
- Changing the iron source controlled cell proliferation rate, which helped maintain appropriate VCD.
- Improvements in the media composition, including L-Cystine, enhanced productivity.

https://www.ajinomotocellistkorea.com/en/technicalmaterial_en



Related Technical Documents (QR-Code)

Customer Service

Global CELLiST™ Solution Centers



CELLiST | Solution Center Korea

70, Songdogwahak-ro, Yeonsu-gu, Incheon, Republic of Korea
Tel: +82-32-210-2695 Fax: +82-32-210-2607
Mail: yaron.silberberg.dk2@asv.ajnomoto.com

CELLiST | Solution Center India

2nd Floor, Block 05, TIMES CORPORATE PARK Thaltej, Ahmedabad, Gujarat 380058, India
Tel: +91-87359-10813
Mail: parshadg.hirapara.ut3@asv.ajnomoto.com

CELLiST | Solution Center China

Room 501, Building 9, NO.185 Jiangnan East Road, Songjiang District, Shanghai, China
Tel: +86-21-6219-2705
Mail: mingqi_fan@ajnomoto.com.cn

For quotations, orders and product samples requests, please contact us at:
CELLiST™ Sales Team
Tel: +82-32-210-2600 Mail: sales@ajnomotocellist.com


Liquid Media Preparation

| CELLiST™ CHO Perfusion Media_P1 (1,000 mL Preparation)

1. Add 37.4 g of CELLiST™ CHO Perfusion Media_P1 to 800 mL cell-culture-grade water (room temperature).
2. Wash the remaining product into the container using cell-culture-grade water.
3. Stir until all powder is dissolved.
4. Add 3.2 g of sodium-bicarbonate.
5. Stir until all particles are completely dissolved.
6. Adjust the volume to 1,000 mL using cell-culture-grade water and stir until the mixture is uniform.
7. Filter the medium in a biosafety cabinet, using a membrane filter with pore size of 0.2 to 0.22 μm in diameter.
(Reference value - pH: 7.0-7.4 / Osmolality: 340-360 mOsm/kg)
8. Store in a refrigerated (2–8°C), dark environment until use.

| L-Cystine (200 mL Preparation)


1. Add 2.1 g of L-Cystine to 160 mL cell-culture-grade water (room temperature).
2. Wash the remaining product into the container using cell-culture-grade water.
3. Add 2.5 mL of 8N NaOH.
4. Stir until all powder is dissolved.
5. Adjust the volume to 200 mL with cell-culture-grade water and stir until the mixture is uniform.
6. Filter the medium in a biosafety cabinet, using a membrane filter with pore size of 0.2 to 0.22 μm in diameter.
(Reference value - pH: 11.5-13.0 / Osmolality: 130-150 mOsm/kg)
7. Store in a refrigerated (2–8°C), dark environment until use.

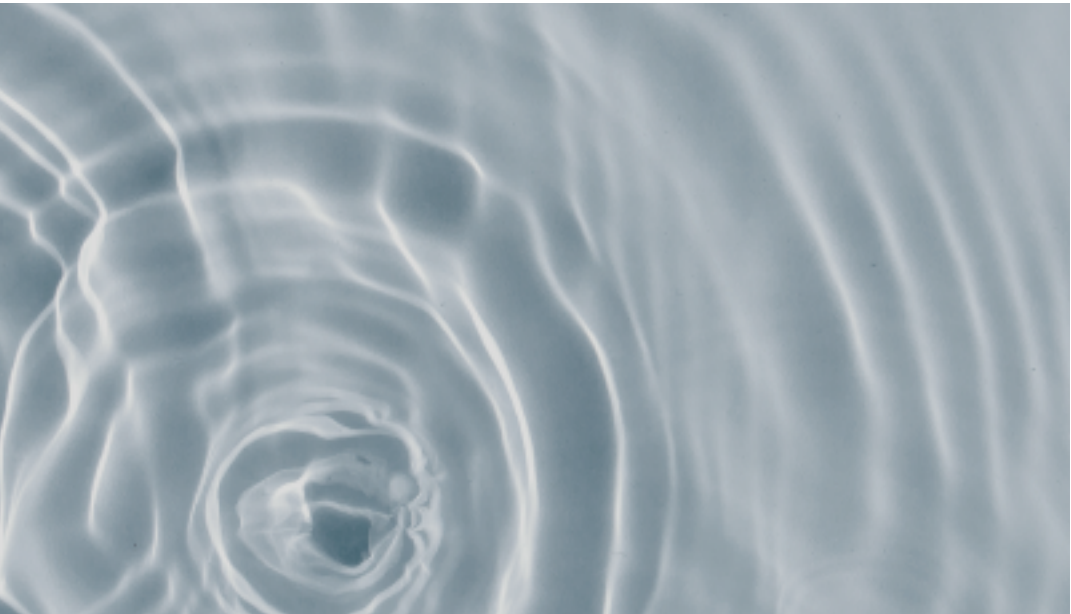

| Storage Conditions

- Before liquid preparation, store powder media in dark, refrigerated environment (2–8°C), away from high humidity.
- After liquid preparation, store in dark, refrigerated environment (2–8°C).

| Usage Guidelines

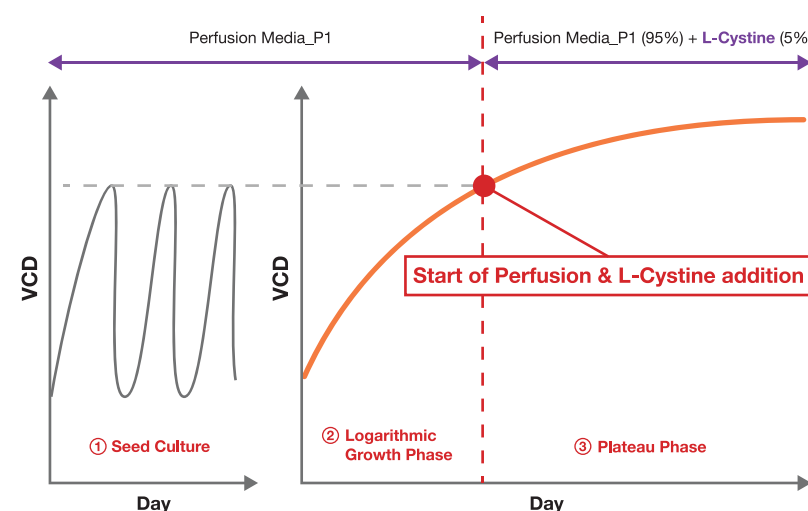

- Right before use, aseptically add 2–6 mM L-glutamine source to CELLiST™ CHO Perfusion Media_P1.
- *Note: CELLiST™ CHO Perfusion Media_P1 can also be used for seed cultures.
- For main culture, use CELLiST™ CHO Perfusion Media_P1 prior to starting perfusion.
- When VCD reaches maximum VCD during seed culture, L-Cystine addition is required, in addition to Perfusion Media_P1, for maintaining adequate culture performance. At the same time, perfusion should also be initiated. We recommend using a ratio of 95:5 (Perfusion Media_P1 : L-Cystine).
- CELLiST™ CHO Perfusion Media_P1 and L-Cystine should not be mixed together directly.

*Note: This product is a cell culture medium for use in research applications. Do not use it for any other purpose.



Eat Well, Live Well.



www.ajinomotocellistkorea.com